EN DIGITAL THERMOMETER model T-30 USER MANUAL

- Measuring time: 20 seconds - oral, 10 seconds - rectal; 20 seconds - axillary
- Last measurement memory
- Waterproof, flexible tip


Thank you for choosing **TM Dr. Frei®** Digital Thermometer Model T-30. We are sure that having apprised worthy the high quality and reliability of this device you will become a regular user of the products of Swiss Trademark **Dr. Frei®**. Before starting to use the thermometer please study the instructions carefully.

IMPORTANT SAFETY INSTRUCTIONS

1. The thermometer is suitable exclusively for measurement of the body temperature!
2. The minimum measurement time until the signaling tone «beep» must be maintained without exception!
3. The thermometer contains small parts (battery, etc.) that can be swallowed by children. Therefore never leave the thermometer unattended to children.
4. Do not store the unit where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures, dropping or subjecting your thermometer to strong shocks.
5. Do not attempt to disassemble the unit, except to replace the battery.
6. Avoid ambient temperatures above 55 °C. NEVER boil the thermometer!
7. Store the unit in the protective case when not in use.
8. Избегайте изгиба наконечника более чем на 45°С.

Batteries and electronic instruments must be disposed of in accordance with locally applicable regulations, nor with domestic waste.

FEATURES OF YOUR DIGITAL THERMOMETER

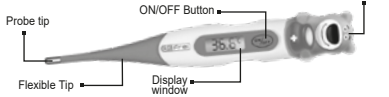HOW TO TAKE YOUR TEMPERATURE

1. Press ON/OFF button to activate. The unit would beep and display **99.9**. This LCD display test will go on for around 2 seconds.
2. Then, the temperature taken last time would be displayed with a «M» mark in the right corner.
3. Then, 37.0°C will appear, when Lo and flashing °C displays, the thermometer is ready for temperature measuring. Thermometer displays Lo under condition the room temperature is lower than 32.0 °C.

4. Put the thermometer to places you want to take temperature, into mouth, under armpit or inside the anus.
5. The beep will sound when temperature measuring is complete (when taking an oral and rectal temperature). When taking an axillary temperature, ignore the beep sound and wait not less than 90 seconds before taking the reading.
6. Press ON/OFF button to turn off the unit. (The unit will automatically power off after approximately eight minutes after use.)

MEASURING METHODS

Measuring method	Position	Measuring time
Taking an oral temperature (in the mouth). Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air. Do not eat or drink cold or hot products before measurement.		20 seconds
Taking a rectal temperature (in the anus). This is the most reliable measuring method, and is especially appropriate for infants and small children. Apply a water-soluble lubricant to the probe cover. Gently insert the probe (MAXIMUM 1/2") into the rectum.		10 seconds
Taking an axillary temperature (under the arm). From a medical viewpoint, this method will always provide inaccurate readings, and should not be used if precise measurements are required.		20 seconds

Note: Rectal temperatures are generally 0.5°C higher than oral readings. Axillary temperatures are generally 0.5°C lower than oral readings.

CLEANING INSTRUCTIONS

Before and after each use, clean the thermometer with a soft cloth and 70% isopropyl alcohol diluted with water, or with COLD soapy water. Do not soak for an extended period of time or sterilize by boiling, gas or steam autoclaves.

REPLACING THE BATTERY

- The battery in your digital thermometer needs to be replaced when a «M» appears in the top right corner of the display. Replace the battery with an AG3 alkaline or LR41 silver oxide type or equivalent. To replace the battery,
1. Pull the battery cover away from the unit.
 2. Use a non-metal pointed tool to remove the battery and insert the new battery

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- with the positive(+)/side facing up.
 3. Carefully slide the battery cover back into place.
- Note:** Please properly replace the battery away from small children and heat.

SPECIFICATIONS

Measurement range	32.0°C ~ 42.9°C
Low temperature display	Temperature > 32.0°C: display: Lo°C
High temperature display	Temperature > 42.9°C: display: Hi°C
Display resolution	0.1°C
Accuracy	35.0°C ~ 39.0°C ± 0.1°C Other range: ± 0.2°C
Fever beep sound	≥ 37.8°C 30 short beep sounds in 10 seconds
Normal beep sound	< 37.8°C 10 long beep sounds in 10 seconds
Auto power off	8 minutes ± 2minutes
Battery	AG3 or LR41
Life of battery	2 years
Operating condition	Temperature: 10°C ~ 40°C Humidity: 30% ~ 85%RH
Storage condition	Temperature: -20°C ~ 55°C Humidity: 25% ~ 95%RH
Water-proof grade	IP22
Dimensions & Weight	140 x 23 x 18 mm, approx. 10 g (including battery)

* The specifications and the outer design are subject to change for improvement reasons without previous notice.

CHECKING

The device is calibrated at manufacturing checking. During usage and after repair it is recommended to check the device at least once a year at a service center in your country.

WARRANTY

The digital thermometer is guaranteed for 24 months from date of purchase. The guarantee does not apply to damage caused by improper handling, accidents, not following the operating instructions or alterations made to the instrument by third parties. The warranty is only valid upon presentation of the warranty card with purchase date and stamp of trade organization.

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