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
# **DEAR DR. FREI MULTIFUNCTIONAL INFRARED THERMOMETER OWNER**

---

Thank you for choosing TM **Dr. Frei®** Multifunctional Infrared Thermometer Model MI-100. We are sure that having appraised worthily the high quality and reliability of this device you will become a regular user of the products of Swiss Trademark **Dr. Frei®**.

This innovative medical device uses advanced infrared (IR) technology to measure temperature instantly and accurately on the body or object.

Before starting to use this device please study the user's manual carefully. The user's manual offers all information you need to correctly measure your blood pressure and pulse. For all questions concerning the device please contact your local distributor or **Dr. Frei®** service center in your country.
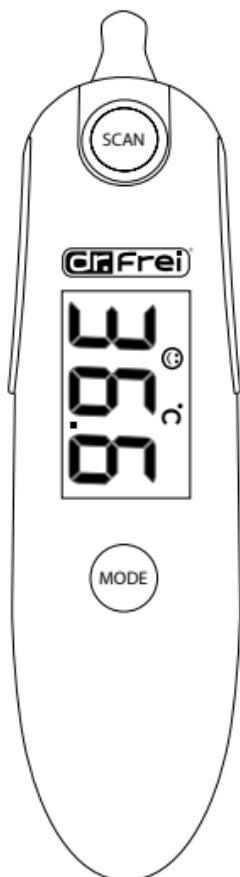
# TABLE OF CONTENTS

---

<b>IMPORTANT INFORMATION ON BODY TEMPERATURE AND ITS MEASUREMENT .....</b>	<b>51</b>
About Normal Body Temperature & Fever .....	51
<b>ADVANTAGES OF MULTIFUNCTIONAL INFRARED THERMOMETER MODEL MI-100 .....</b>	<b>52</b>
Important information .....	52
<b>BEFORE YOU BEGIN .....</b>	<b>53</b>
Cautions and Warnings .....	65
Meter Overview and Display Screen .....	53
Setting the time.....	55
Selecting the Temperature Unit .....	55
Select the Display Mode.....	55
Replacing the Battery .....	57
<b>METER MEMORY .....</b>	<b>58</b>
Recalling the Memory.....	58
<b>ERROR MESSAGES .....</b>	<b>58</b>
<b>MAINTENANCE .....</b>	<b>59</b>
<b>SYMBOL INFORMATION .....</b>	<b>60</b>
<b>SPECIFICATIONS.....</b>	<b>61</b>
<b>WARRANTY .....</b>	<b>62</b>

# **IMPORTANT INFORMATION ON BODY TEMPERATURE AND ITS MEASUREMENT**

---

## **About Normal Body Temperature & Fever**

---

The temperature in the forehead and temple area differs from the internal temperature, which is taken orally or rectally. Vasoconstriction, an effect which constricts the blood vessels and cools the skin, can occur during the early stages of a fever. If the measurement therefore does not match the patient's own perception or is unusually low, repeat the measurement every 15 minutes. As a reference, you can also measure the internal body temperature using a conventional oral or rectal thermometer.

Body temperature can vary from one individual/person to next. It also varies by location on the body and time of day. Below shows the statistical normal ranges from different sites. Please keep in mind that temperatures measured from different sites, even at the same time, should not be directly compared. Fever indicates that the body temperature is higher than normal. This symptom may be caused by infection, overdressing or immunization. Some people may not experience fever even when they are ill. These include, but are not limited to, infants younger than 3 months old, persons with compromised immune systems, persons taking antibiotics, steroids or antipyretics (aspirin, ibuprofen, acetaminophen), or persons with certain chronic illnesses. Please consult your physician when you feel ill even if you do not have fever.

**Table \*1 Body Site Normal Temperature Range**

<b>Body Site</b>	<b>Normal Temperature Range</b>
Oral	35.5°C-37.5°C (95.9°F - 99.5°F)
Ear	35.8°C-38°C (96.44°F - 100.4°F)
Axillary	34.7°C - 37.3°C (94.46°F - 99.14°F)

\* Body Temperature at WebMD;  
website: <http://firstaid.webmd.com/body-temperature> retrieved at 2010 Jan 7.

# **ADVANTAGES OF MULTIFUNCTIONAL INFRARED THERMOMETER MODEL MI -100**

---

## **Important information**

---

This thermometer is clinically proven to produce accurate temperature measurements.

The thermometer model MI-100 is intended for the intermittent measurement and monitoring of human body temperature from forehead or ear measurement at home.

The thermometer measures the infrared heat generated by the surface of the skin over the vessel and its surrounding tissue.

The thermometer then converts it into a temperature value.

- Auto standby if not in use for 25 seconds will auto skip to time display.
- Memory function allows you to store 10 most recent body mode temperature results.
- Easy to read LCD with backlight in a dark environment.
- Color visible indication of alarm slight fever (red) and normal reading (green).

### **NOTE:**

Please be advised that the accuracy cannot be ensured when the thermometer is not clean. Check that the detecting sensor is clean before taking a measurement

A control measurement using a conventional thermometer is recommended in the following cases:

1. If the reading is surprisingly low.
2. For new-born infants up to 100 days old.
3. For children under three years of age who have a weakened immune system or who react unusually in the presence or absence of fever.

The thermometer does not emit any infrared signal.

# BEFORE YOU BEGIN

---

## Cautions and Warnings

---

1. As with any thermometer, proper technique is crucial to obtaining accurate temperature readings. Please read this manual thoroughly and carefully before using.
2. Always operate the thermometer in an operating temperature range 15°C to 25°C (59°F to 77°F), and relative humidity less than 85%.
3. Always store the thermometer in a cool and dry place -25°C to 60°C (-13°F to 140°F) and relative humidity less than 85%, atmospheric pressure: 70kPa~106kPa.
4. Avoid direct sunlight.
5. Avoid dropping the thermometer, if it happens and you think the thermometer may be damaged, please contact customer services immediately.
6. Do not disassemble the thermometer.
7. Please don't touch or blow the infrared sensor. Damaged or soiled infrared sensor may lead to incorrect measurement results.
8. Please do not insert infrared sensor into the ear forcibly, please stop using upon you have the abnormal feelings such as pain.  
\*Or it could hurt the ear canal.
9. The people who are suffering from otitis externa, otitis media and so on ear diseases, please do not use the ear mode.  
\*Otherwise it may cause the deterioration of injured part.
10. When the ear canal is wet after swimming or shower. Please do not use.  
\*Or you may hurt the ear canal.
11. Please don't use this product in any other body part except ear and forehead.

## Meter Overview and Display Screen

---

### Package Contents:



Thermometer

User's Manual

2\*1.5V AAA battery

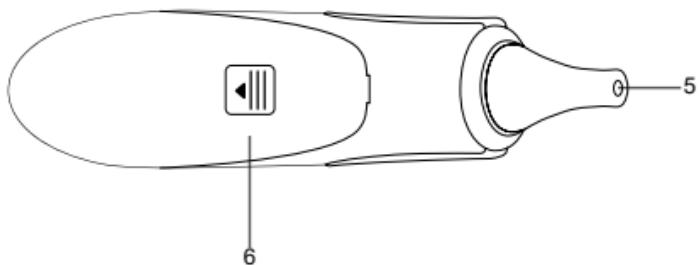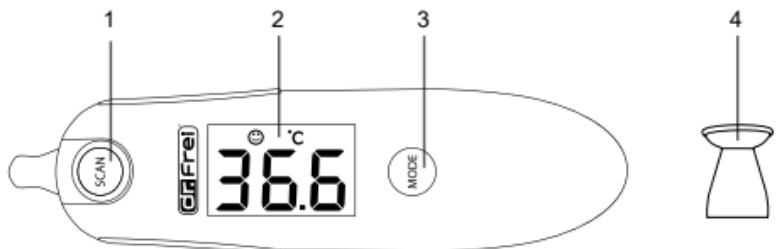1. On/Scan Button  
 2. Display Screen  
 3. Mode Light  
 4. Forehead temperature measuring attachment  
 5. Detecting Sensor  
 6. Battery Compartment

1. Environment temperature indicator  
 2. Forehead temperature indicator  
 3. Ear temperature indicator  
 4. Clock  
 5. Face symbol  
 6. Temperature unit  
 7. Memory  
 8. Temperature reading

## **Setting the time**

---

Take out the batteries and then replace them. The display will show "12:00" and "12" will flash after two seconds. Adjust the hour pressing SCAN button (1) and fixing the correct option using MODE button (3).

"00" will flash on the display after you fix the hour. Adjust minutes pressing SCAN button (1) and fixing the correct option using MODE button (3).

## **Selecting the Temperature Unit**

---

When device is ON, press SCAN button until the display screen shows the correct temperature symbol for your use, exactly °C or °F. Then press the SCAN button once again to complete the setting.

## **Select the Display Mode**

---

### **1. Body Mode.**

Ear temperature measurement

- When device is ON, press the SCAN button (1). After two beeps the last measurement will be displayed for 3 seconds together with the °C or °F flashing. The icon is shown.
- Insert infrared sensor into the ear and press SCAN button (1). The sensor probe will automatically scan the ear membrane's temperature after a beep.
- When measurement is complete, the result will be shown on the display screen after one long beep.
  - \* Insert the sensor probe into the ear canal as deeply as possible, WITHOUT causing pain.
  - \*\* If your ear canal is too small and it is difficult to insert the probe, pull back the ear gently and insert infrared probe into the ear canal carefully and WITHOUT any force.

Forehead Temperature Measurement

- Attach the forehead temperature measuring attachment (3) to the sensor probe (1) of the thermometer before use.

- b) When device is ON, press the SCAN button (1). After two beeps the last measurement will be displayed for 3 seconds together with the °C or °F flashing.
- c) Press the MODE button (3) to switch the thermometer into forehead measuring mode (icon).
- d) Press SCAN button (1). The sensor probe will automatically scan the forehead temperature after a beep.
- e) When measurement is complete, the result will be shown on the display screen after one long beep.

#### Additional Indicators

Ear and forehead temperature measurements:  $\leq 37.4^{\circ}\text{C}$  - the display screen shows smile face symbol,  $\geq 37.5^{\circ}\text{C}$  - the display screen shows crying face symbol Backlight indicator.  $\leq 37.4^{\circ}\text{C}$  ( $99.32^{\circ}\text{F}$ ) - the display screen turns GREEN,  $\geq 37.5\sim 37.9^{\circ}\text{C}$  ( $99.5^{\circ}\text{F} \sim 100.22^{\circ}\text{F}$ ) - the display screen turns ORANGE,  $\geq 38.0^{\circ}\text{C}$  ( $100.4^{\circ}\text{F}$ ) - the display screen turns RED

#### 2. Object or Ambient Temperature Mode.

##### Object Temperature Measurement

- a) When device is ON, press the SCAN button (1). After two beeps the last measurement will be displayed for 3 seconds together with the °C or °F flashing.
- b) Press the MODE button (3) twice to switch the thermometer into object temperature measuring mode (without any icon displayed).
- c) Press the SCAN button (1) again keeping 10 mm distance between the thermometer and the object measured.
- d) The display screen will display “---” after one beep and make temperature sweep of the object measured automatically.
- e) The result of measuring will be displayed on the screen after one long beep.

##### Ambient Temperature Measurement


- a) When device is ON, press the SCAN button (1) to turn on. After two beeps the last measurement will be displayed for 3 seconds together with the °C or °F flashing.
- b) Press the MODE button (3) three times to switch the thermometer into ambient temperature measuring mode (icon is displayed).


## Replacing the Battery

The thermometer comes with two 1.5V alkalescence batteries (AAA). The meter will display «» to alert you when the meter power is getting low, please follow the steps below to replace new battery immediately.

1. Press the buckle on battery cover and lift up to remove the cover.
2. Remove the old battery and replace with two 1.5V alkalescence batteries (AAA).

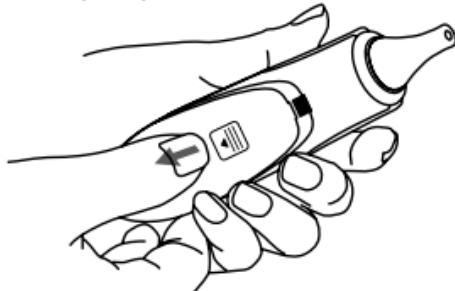3. Close the battery cover.
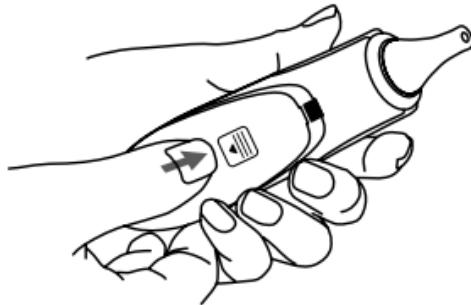
### NOTE:

1. Although the thermometer works when «» appears, we still recommend that you change the battery to obtain an accurate result.
2. Remove the battery if stored for a long period of time.
3. The battery should be kept away out of child's reach.
4. The battery shall be properly disposed according to your local environmental regulations.
5. Follow the national requirement to dispose the unit.

# METER MEMORY

## Recalling the Memory

Your thermometer stores 10 most recent body temperature readings.

Previous Measurement Temperature

- When device is ON, press the SCAN button (1) to turn on.
- Keep the MODE button (3) pressed to enter the list of previous measurements for 1 second, "M" will be displayed.
- Every time you press the MODE button (3), a result will be displayed in the order.

## ERROR MESSAGES

If you follow the recommended action but the problem persists, or error messages other than the ones below appear, please call your local distributor or Dr. Frei® service centre in your country. Do not attempt to repair by yourself and never try to disassemble the monitor under any circumstances.

Error display	Reason	Solution
No display	Batteries incorrectly inserted (+ve/-ve).	Install the batteries correctly
	Batteries empty	Replace with new batteries

Confused result of measurement	Measurement not continuous	Measure temperature after 3 minutes
	Sensing head is dirty	Use a cotton swab or soft cloth and wipe it gently with alcohol
	Temperature sensing head is damaged or scratched	Please contact agent to repair. (Must be under warranty)
	The ear has foreign matter within it.	Consult a doctor
	Measurement done incorrectly	Please read the manual carefully
Display shows "HH"	Indicates the system detects the ambient temperature exceeds the normal operational range of the device	The measurement should be done in an environment between 15°C to 20°C and the device should be stabilized for 30 minutes in such an environment
Display shows "LI"	Indicates the system detects the ambient temperature exceeds the normal operational range of the device	The measurement should be done in an environment between 15°C to 20°C and the device should be stabilized for 30 minutes in such an environment
Display shows "Hi"	Indicates that the target temperature is above the measurement limit of the device	
Display shows "Lo"	Indicates that the target temperature is above the measurement limit of the device	

## MAINTENANCE

---

### Cleaning and Care

---

- The detecting sensor is not waterproof. Please wipe with a

- clean and dry cotton swab to clean the infrared sensor.
- The body of the thermometer is not water-resistant. Never put the thermometer under a running tap or submerge it into water. Use a soft and dry cloth to clean it. Do not use abrasive cleaners.
  - Store the thermometer in a cool and dry location. Free from dust and away from direct sunlight.

## SYMBOL INFORMATION

---

Symbol	Referent
	Consult instructions for use
	Manufacturer
<b>SN</b>	Serial number
	Caution, consult accompanying documents
<b>CE0123</b>	Mark CE
	Authorized representative in the European Community
	Type B Equipment
	Collection for electrical and Electronic equipment

## SPECIFICATIONS

<b>Model No</b>	MI-100
<b>Dimensions &amp; Weight</b>	145 x 37 x 35 mm, 53 g (include batteries)
<b>Power Source &amp; Voltage</b>	dc 3 V (Two 1.5V alkalescence batteries (AAA))
<b>Displayed temperature range</b>	0°C to 100°C (32°F to 212°F)
<b>Display resolution</b>	± 0.1°C/0.1 °F
<b>Accuracy</b>	<b>Ear mode:</b> ±0.2°C (35.0°C - 42.0°C)/ ±0.4°F (95.0°F ~ 107.0°F)  <b>Forehead mode:</b> ± 0.3°C (34.0°C ~ 44.0°C)/ ± 0.5°F (<95.0°F / > 107.6°F)
<b>Temperature unit</b>	°C or °F
<b>Operating temperature range</b>	15°C to 25°C (59°F to 77°F)
<b>Operating humidity</b>	30% RH - 85% RH
<b>Storage temperature range</b>	-25°C to 60°C (-13°F ~ 140°F)
<b>Storage humidity</b>	10% RH - 85% RH
<b>Memory capacity</b>	10 measurements
<b>Standards</b>	Device Corresponds to the requirements of the standard for infrared thermometers ASTM E1965-98  EN 12470-5 EN 60601-1-2:2007/AC:2010 EN 60601-1-4:1996 EN 60601-1-6:2010

	<p><b>Electromagnetic Compatibility:</b> Device fulfills the stipulations of the standard EN 60601-1-2</p> <p>The stipulations of EU-Directive 93/42/EEC for Medical Devices Class IIa have been fulfilled.</p>
--	---

\* The specifications and the outer design are subject to change for improvement reasons without previous notice.

Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, Not with domestic waste.

## **WARRANTY**

---

This product is under warranty for 2 years after date of purchase. The warranty is only valid upon presentation of the warranty card which was correctly filled in and sealed by the seller (or by the service center representative) with the date of the purchase (repair), the check, instruction manual and safe packaging of the device.

- The warranty does not apply to wear out of battery and packaging box.
- The warranty does not apply to damage caused by improper handling, accidents, not following the operating instructions or self-maintained alterations made to the device.
- The device cannot be returned if it is without original packaging, with mechanical damages (scratches, stains and so on) or without any component.
- The device which was used and needs repair cannot be replaced with a new one.

In case there was violation of operating rules during warranty period, repair will be done at the expense of the consumer.

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